



CDP SWISS
MEDICAL
PRODUKTKATALOG

inspiring customers worldwide

ROBOLEX




ROBOLEX

ROBOLEX is an effective and space saving body contouring device combining Vacuum suction, Ultrasound Cavitation, Low Level Laser and Multi/Bi-Polar Radio Frequencies.

Lipocyte destruction by using Ultrasound Cavitation.

Dissolving destroyed lipocyte by using Multi/Bi-Polar Radio Frequency and cellulite reduction by Low Level Laser.


Drain the dissolved lipocyte rapidly from the body to outside by the lymph drainage and blood circulation and maintain the elasticity on treated area with Vacuum.

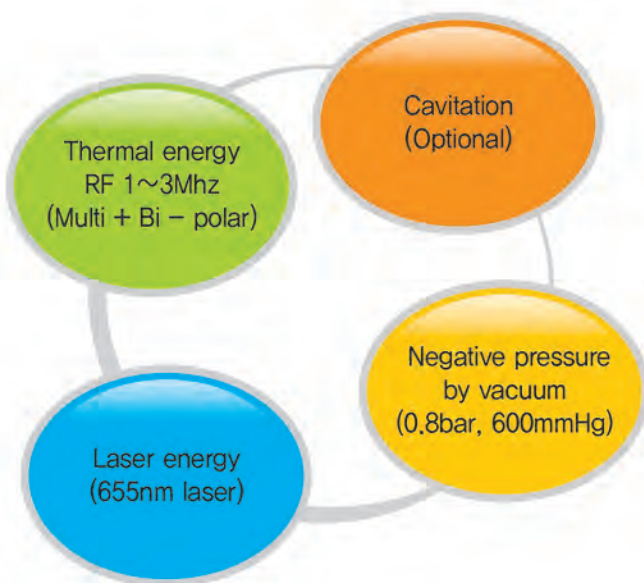

Face/Hair



Arm/Leg/Thigh

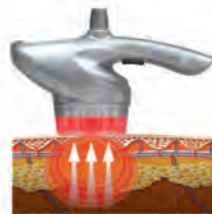

Belly/hip



Body



RV-25



RV-60




RV-90



CAVI-38(Optional)

Dissolve fat where under 5mm from the skin by 655nm laser
 RF infiltrates into 5~15mm from the skin and spread oxygen into intracellular by heating
 Soften skin with negative pressure made by vacuum and help to distribute energy evenly

APPLICATION

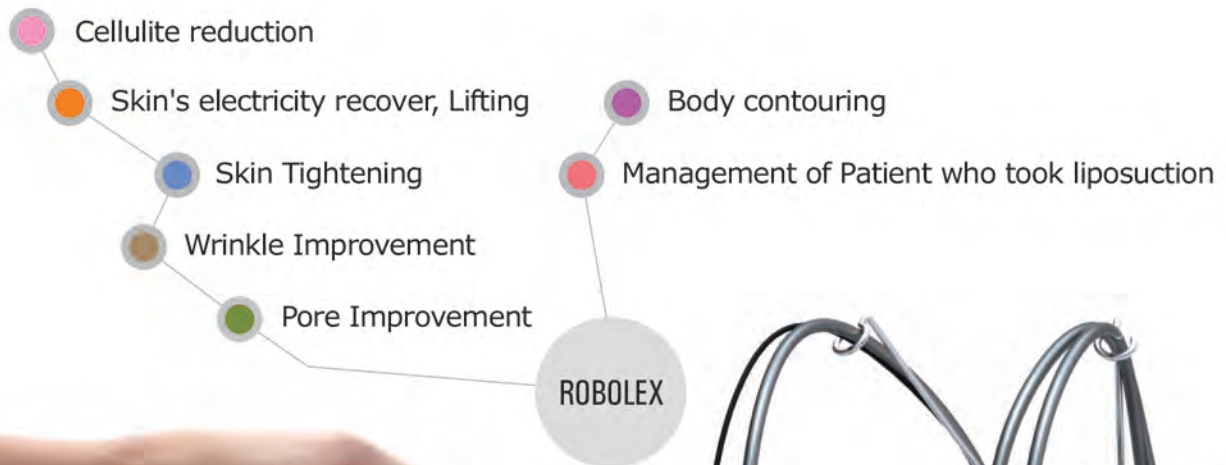

SUMMARY

METHOD

Obesity care by Multi-functional treatment (RF + Endomologie + Laser)
Optional item : Ultrasound cavitation.

FEATURE

Completion of wanted body shape by selective reduction of fat cell and Cellulite with non invasive treatment.

THEORY

Destruction of lipocyte by combination operation with negative pressure in lipocyte & ultrasound cavitation.

SAFETY

Complementary lipolysis treatment way in safe to prevent the risk by liposuction and to get rid of uncertainty by injection therapy.

EFFECT

Confirmation of reduction for size & body shape by 3~4 times treatment and it prevents the wrinkled line on the reduced part after treatment.



ROBOLEX



POWER	AC110 ~ 220V, 50/60Hz
POWER CONSUMPTION	MAX 500VA
RF OUPUT(MAX)	20W(at 500Ω)
RF OUTPUT FREQUENCY	1MHz ~ 1.8MHz(RV95, 60) 2MHz ~ 3MHz(RV25) (Automatic frequency adjustment)
RF OUTPUT TYPE	MULTI(RV95, 60), BI(RV25)
VACUUM	MAX 0.8bar(600mmHg) (AC Pump Cont. Max 2.1Bar)
LASER	655nm(3mW)-RV95,60
ULTRASOUND FREQUENCY (Optional)	38Khz

CDP SWISS^{CH} MEDICAL

PRODUKTKATALOG

Hauptsitz

CDP SWISS AG
Weinbergstrasse 22
CH-8001 Zürich

Niederlassung

CDP SWISS AG
Moosstrasse 31
CH-8907 Wettswil am Albis

Tel. +41 43 538 57 45
Fax +41 43 538 57 73

www.cdpswiss.com